

FineCross® M3

Coronary Micro-Guide Catheter


Runthrough® NS Izanai™

Coronary Guidewires

**UNIQUELY
ENGINEERED**
FOR SIMPLE TO
COMPLEX PCI



 **TERUMO**
INTERVENTIONAL
SYSTEMS

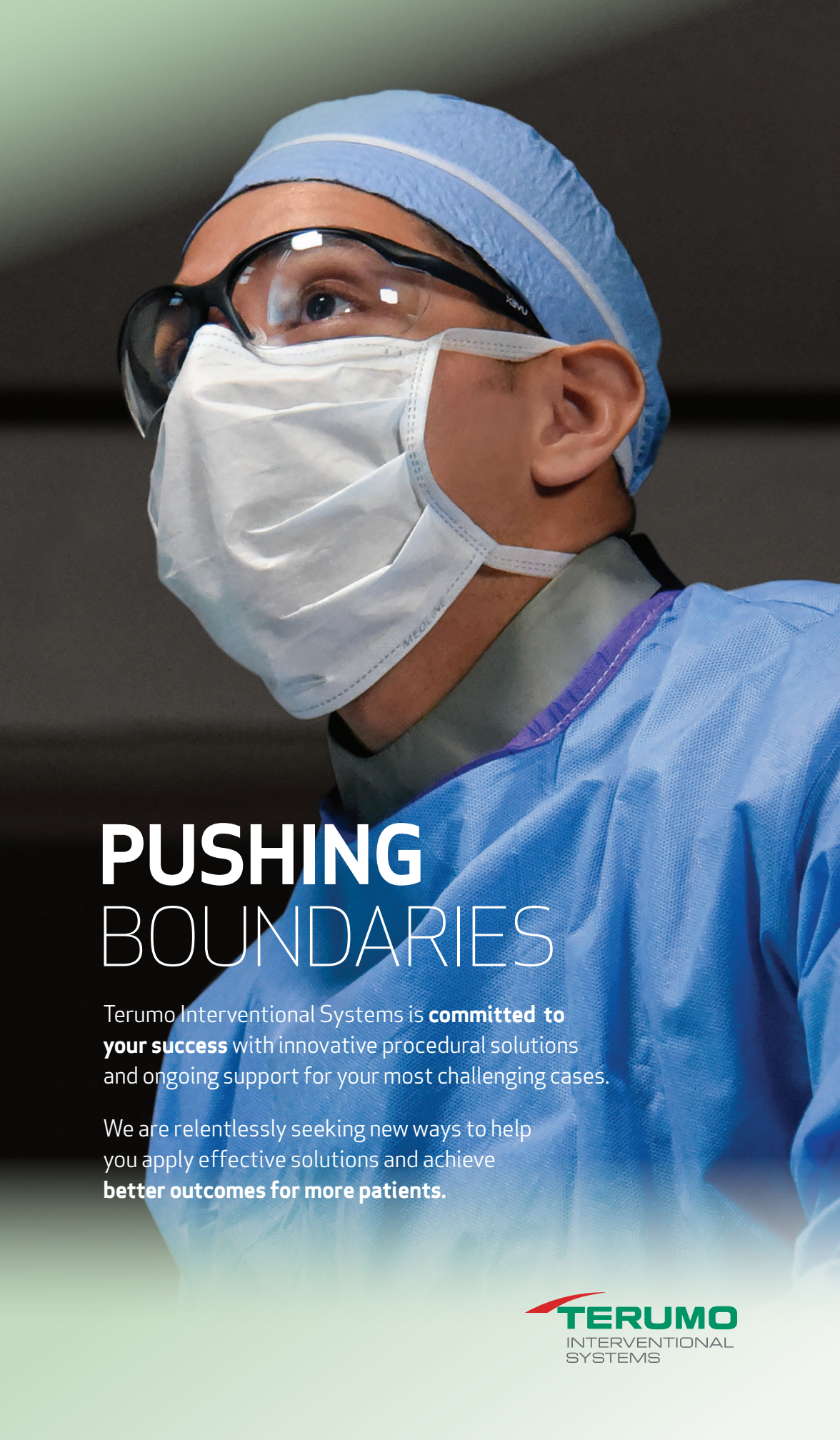

PUSHING BOUNDARIES

Terumo Interventional Systems is **committed to your success** with innovative procedural solutions and ongoing support for your most challenging cases.

We are relentlessly seeking new ways to help you apply effective solutions and achieve **better outcomes for more patients.**




FineCross® M3

Coronary Micro-Guide Catheter

OPTIMAL GUIDEWIRE SUPPORT AND SUPERIOR CROSSABILITY¹

FINECROSS® M3 Coronary Micro-Guide Catheter is designed to provide the ultimate balance of properties for access and crossing of complex lesions.


- 1.7 Fr tapered distal tip is designed for improved crossability and guidewire handling¹
- Next Generation Hydrophilic Coating, with its unmatched and most durable lubricity, provides excellent trackability and crossability¹
- 15 cm floppy distal segment allows for improved trackability
- Tapered stainless steel braid construction is designed for optimal guidewire support

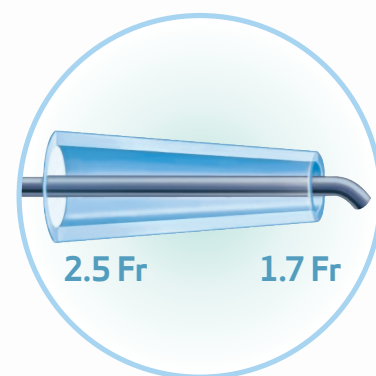

Tapered Distal Tip

Tapered ID and OD along the length of the catheter for unsurpassed crossability²

Outer proximal diameter: 2.5 Fr
Outer distal diameter: 1.7 Fr

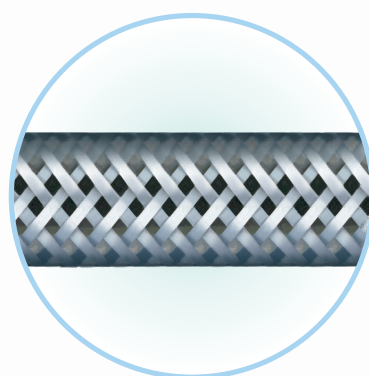

Next Generation Hydrophilic Coating

Durable lubricity provides excellent crossability and trackability



Stainless Steel Braided Mesh

Provides guidewire support for improved pushability

Floppy Distal Segment

Distal 15 cm flexible length improves trackability through tortuous anatomy

PTFE Inner Layer

Reduces friction and improves guidewire handling

Gold Radiopaque Marker

For enhanced visibility that facilitates accurate placement

Soft, Atraumatic Distal Tip

1.7 Fr entry profile, tapered distal tip for optimal guidewire handling and improved crossability

PROVEN TECHNIQUES TO OPTIMIZE YOUR APPROACH

Access and cross simple to complex lesions using FINECROSS® M3 Coronary Micro-Guide Catheter with RUNTHROUGH® NS Coronary Guidewires


SUPPORT for Side-branch or Acute Bifurcation Lesions

Use **FineCross® M3**

Coronary Micro-Guide Catheter

with **Runthrough® NS Extra Floppy**

Coronary Guidewires

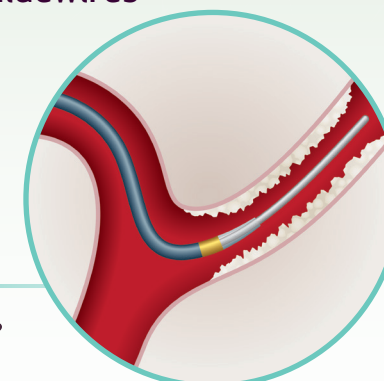


TRACKABILITY for Distal or Torturous Lesions

Use **FineCross® M3**

Coronary Micro-Guide Catheter

with **Runthrough® NS Izanai™**

Coronary Guidewires

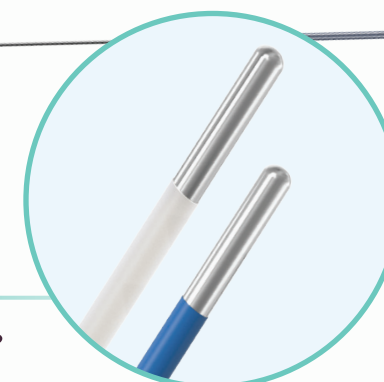

RUNTHROUGH NS Izanai is available in blue and white for simple identification of secondary wires


CROSSABILITY for Tight Stenotic or Calcified Lesions

Use **FineCross® M3**

Coronary Micro-Guide Catheter

with **Runthrough® NS Hypercoat™**

Coronary Guidewires

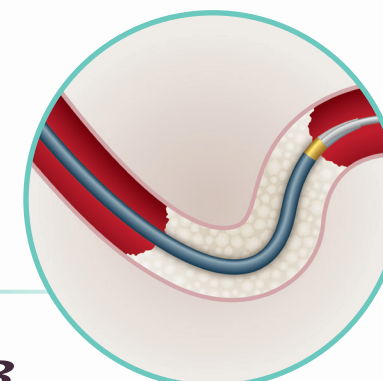

TERUMO
INTERVENTIONAL
SYSTEMS

FineCross® M3

Coronary Micro-Guide Catheter

Runthrough® NS

Coronary Guidewires

ORDERING INFORMATION

FINECROSS® M3 Coronary Micro-Guide Catheter

PRODUCT CODE	EFFECTIVE CATHETER LENGTH (cm)	DISTAL ENTRY TIP PROFILE (FR/MM)	DISTAL OUTER DIAMETER (Fr/mm)	PROXIMAL OUTER DIAMETER (Fr/mm)	DISTAL HYDROPHILIC COATING (CM)	DISTAL INNER DIAMETER (MM)	PROXIMAL INNER DIAMETER (IN/MM)	MAXIMUM GUIDEWIRE SIZE (IN/MM)
M3-1430	130	1.7/0.55	1.8/0.61	2.5/0.84	75	0.45	0.021/0.55	0.014/0.36
M3-1450	150	1.7/0.55	1.8/0.61	2.5/0.84	95	0.45	0.021/0.55	0.014/0.36

RUNTHROUGH® NS Extra Floppy Coronary Guidewire

PRODUCT CODE	DIAMETER (in)	LENGTH (cm)	TIP RADIOPACITY (cm)	TIP SHAPE
25-1011	0.014	180	3	Straight/Shapeable
25-1013	0.014	300	3	Straight/Shapeable

RUNTHROUGH® NS Izanai™ Coronary Guidewire

PRODUCT CODE	DIAMETER (in)	LENGTH (cm)	TIP RADIOPACITY (cm)	TIP SHAPE / SHAFT COLOR
25-5011	0.014	180	3	Straight/Blue
25-5211	0.014	180	3	Straight/White
25-5213	0.014	300	3	Straight/White

RUNTHROUGH® NS Hypercoat™ Coronary Guidewire

PRODUCT CODE	DIAMETER (in)	LENGTH (cm)	TIP RADIOPACITY (cm)	TIP SHAPE
25-3011	0.014	180	3	Straight/Shapeable
25-3013	0.014	300	3	Straight/Shapeable

Packed 5 per box. Shaping tool included.

RUNTHROUGH® NS Extension Wire

PRODUCT CODE	DIAMETER (in)	LENGTH (cm)	TIP RADIOPACITY (cm)	TIP SHAPE
25-9092	0.014	150	—	—

Packed 1 per box. Does not include a shaping tool.

FIND OUT MORE  Phone: 800-888-3786  [terumo.com](https://www.terumo.com)

References:

1. Data on file.
2. Claims limited to tested microcatheters. (Volcano Valet & ASAHI Corsair)

RX ONLY. Refer to the product labels and package insert for complete warnings, precautions, potential complications, and instructions for use.

©2023 Terumo Medical Corporation. All rights reserved. All brand names are trademarks or registered trademarks of their respective owners. PM-06949

